CURAPROX

HYDROSONIC
EASY

Manual

 SWISS PREMIUM ORAL CARE

CURAPROX

HYDROSONIC
EASY

Deutsch

ab Seite 6

English

see page 12

Français

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Italiano

vedi pagina 24

Español

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Português

a partir da página 36

Česky

od strany 42

Slovenščina

Pojdite na stran 48

Slovensko

Glejte stran 54

Warnhinweis

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HOW TO USE

Warning

If you've recently had surgery on your teeth or gums, consult your dentist before using the 'Hydrosonic easy'. If you experience excessive bleeding after using this device, please contact your dentist immediately. This toothbrush can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision. Children shall not play with the appliance. The 'Hydrosonic easy' complies with safety standards for electromagnetic devices. If you wear an implant or pacemaker and have questions about the device, please contact your doctor or the manufacturer of your implant or pacemaker.

Safety note

Warning: For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance. This appliance contains batteries that are non-replaceable. Before using, please check that the voltage stated on the device complies with your local mains voltage. If the Hydrosonic charger or connecting cable is damaged, please do not use or attempt to repair it yourself. The power cord cannot be replaced. If the cord is damaged, the entire device must be replaced. All repairs must be handled by a

dealer or our customer service. Keep the charger away from water, and do not immerse in water or other liquids. Cleaning the charger: Disconnect from the mains before cleaning and reconnect only once the charger is completely dry. Never charge the device unattended.

Environmentally friendly disposal

Do not dispose of this product with normal household waste at the end of its service life. Take it to a recycling centre or official collection point. This helps protect the environment. The integrated rechargeable battery contains substances that may endanger the environment. So please take the used handset to an official recycling centre or collection point.

Proper disposal of this product

Within the EU, this symbol indicates that this product must not be disposed of with household waste. Old appliances contain valuable recyclable materials. They should be recycled in

order to prevent harm to the environment or human health through improper waste disposal. Therefore, please dispose of old devices through a proper collection system, or send the device back to where you purchased it for proper disposal. They will then ensure that the device is recycled accordingly.

Technical data for USB plug

Adapter input: 100-240 VAC, 50/60Hz, 0.2A max.

Adapter output: 5.0 VDC, 600mA

Technical data for charging station with USB cable

Input: 5.0 VDC, 600mA


Technical data for handset

- Battery: lithium-ion (3.7 VDC, 700mAh)
- Charge the 'Hydrosonic easy' completely before first use. Optimal battery charging time: 10 hours
- Operating time per battery charge: approx. 60 minutes. This is equivalent to 15 days if you brush your teeth for four minutes each day.
- Brushing movements per minute: mode «start» 44'000, mode «clean» 64'000, mode «smile» 84'000.
- Quadrant timer: 30 seconds
- Brushing timer: two minutes
- Automatic switch-off: after four minutes

Contents

- 1× 'Hydrosonic easy'
- 1× 'sensitive' brush head
- 1× charger with USB cable
- 1× USB plug
- 1× travel case
- 1× instruction manual

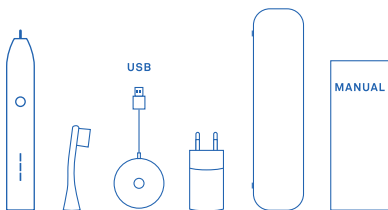
QUICK START

1. Attach brush head

Caution: Do not bend.

2. Switch on

Press button briefly.

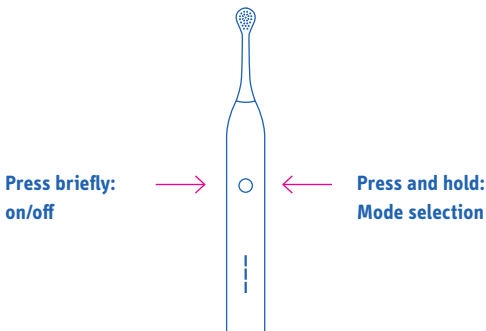3. Mode selection

Press the button and hold it down to switch to the next mode: from «start» to «clean» to «smile». When you reach the desired mode, release the button.

4. Brushing your teeth

Clean from tooth to tooth, fully relaxed; let the bristles do the work; keep the brush as still as possible, pressing down only slightly, or ideally not at all.

5. Two-minute timer with 30-second interval

- After 30 seconds, three short, consecutive pauses of 0.2 seconds each will indicate that it's time to move to a different quadrant in your mouth.
- A long interruption indicates that two minutes have passed. After four minutes, the 'Hydrosonic easy' turns off automatically.

6. Switch off

Press button briefly, and the mode you used last will be saved for next time.

CUSTOMER SERVICE

Guarantee terms

Curaden AG provides guarantee for material, function and manufacturing defects for two years from the date of purchase. This may be carried out by way of repair work or a replacement device. In order for a guarantee claim to be made, the original purchase receipt must be submitted. Keep your receipt as proof of purchase. The guarantee shall expire if the device is subject to interference or changes to its original state.

Guarantee exclusions

In the following cases, Curaden AG shall not provide guarantee services: Wear and tear of the device or parts of the device due to normal use, or if the device was not used, not properly maintained, or if the safety instructions provided with the device were not observed.

Repairs may only be carried out by service partners authorised by Curaden AG.

Service

In order to make a guarantee claim, send the device and original receipt to your country's official Curaprox representative. If you do not know who your country's representative is, visit our website www.curaprox.com or contact us directly in Switzerland: ask@curaprox.com.

- Please discuss any questions about dental or orthodontic care with your dentist. You can also contact us, especially if you have technical problems or if something is not working: ask@curaprox.com.
- Our FAQs are often helpful: www.curaprox.com/faq. Here, you will also find out how to exchange your 'Hydrosonic easy' if necessary.
- The 'Hydrosonic easy' isn't for you? You can return it for a full refund within 30 days of the purchase date.

HOW TO BRUSH YOUR TEETH WITH YOUR 'HYDROSONIC EASY':

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